


*The Sling Specialist
In Safe Patient
mobilization*



SLING CATALOG
Reusable, Wipeable & SPU Series
July 2018

800.273.5749
www.alphamodalities.com
info@alphamodalities.com

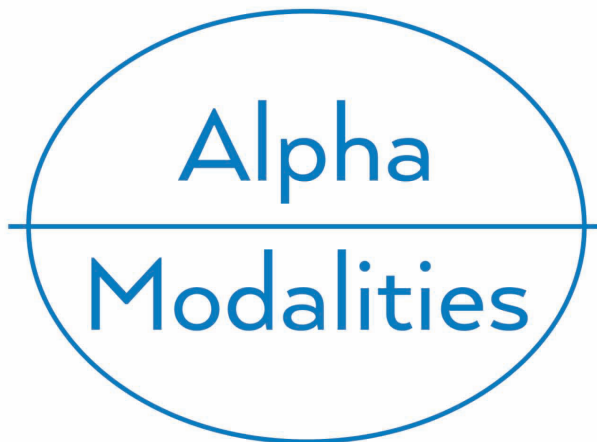

Guide to Color Coding for Mobility Assessment

RED Full Assist: *Patient is unwilling or unable to cooperate*

Page 1.....	A-TSL Series
Page 2.....	A-Seat Series
Page 3.....	A-SitOnSling Series
Page 4.....	A-TTurner Series
Page 5.....	A-Pannus/A-PannusSW
Page 6.....	A-Limb

YELLOW Partial Assist: *Patient can follow basic commands but requires assistance*

Page 7.....	A-SitStandSW
Page 8.....	A-WalkVSW
Page 9.....	A-BariV/A-BariVSW

GREEN Stand-by Assist: *Patient can follow basic commands and requires little to no assistance*

Page 10.....	A-WedgeSW
Page 11.....	A-BedLadderSW

Page 12.....	SPU Series
Page 13.....	Customizations
Page 14.....	Compatibility Statement
Page 15.....	Cleaning and Warranty

Alpha Modalities

A-TSL Series

Intended use: Turns, boosting up in bed, lateral transfers, and cardiac chair transfers.



Product Features

- Breathable—able to remain under patient
- Attaches to ceiling and floor lifts w/ 2,4,6 or 8 point hanger bar systems
- Pressure Mapping documentation on low air loss support surfaces
- Available in multiple colors
- SWL of 454kg /1000lb w/option for 500kg/ 1100lbs
- Latex Free

Options

Materials

- Open Weave Mesh—bathing and general use
- Tight Weave Mesh—replaces linen on bed
- IT—wicking material w/woven antimicrobial

Sizes

- Width Options: 32", 45", 50", 60", or 66"
- Length Options: 64", 86", 100", or 118"

Attachments

- Color Ladder (CL)—industry standard
- F8 Numerical— numerical identification aids caregivers w/color blindness
- F8 Recoil (RC)- designed to prevent tripping with low support surfaces



Alpha Modalities

A-Seat Series

Intended use: Seated transfers to and from bed, chair, toilet and wheelchair with head support.



Product Features

- Tight-weave, open weave bathing mesh and wipeable series.
- Our Bari Line is made with our IT wicking material and will address patients of size up to 500kg /1100 lbs
- We also are introducing a seated sling without binding as pressure point so that it can be left behind the patient for longer periods.

Options

- Available with or without head support
- Loop or key clip attachment
- Color coded to easily identify size

Sizes/Recommended Weight Ranges:

- **Pediatric** (Dino Print): < 15kg
- **Small** (Red): 14-34kg
- **Medium** (Yellow): 34-66kg
- **Large** (Green): 66-100kg
- **X-Large** (Purple): 100-200kg
- **XX Large** (White): 200-363kg
- **Bari-XX** (Grey w purple trim): 363- 500kg



A-SitOnSling Series

Intended use: Designed to remain underneath the patient and provide full body support w/out splitting the legs when transferred. Toileting option available.



Product Features

- Made with our Grey IT material designed to remain underneath the patient for longer periods.
- The edge of sling color is color coded to indicate size
- The sling offers support from the knees to the shoulders on patients.
- Placement is accomplished by turning the patient side to side, positioning the sling at the bend of the knee.
- Slings are washable.
- Latex Free

Options

Sizes

- **Medium** (Yellow Trim) - Seat width 30in/72cm
- **Large** (Green Trim) - Seat width 32 in/81cm
- **XL** (Purple Trim) - Seat width 35in/89cm
- **XXL** (White Trim) - Seat width 39in/99cm

Materials

- Grey IT material
- Polyester Webbing and Trim
- Coated Polyester Thread
- Max SWL 454kg/1000lbs



Alpha Modalities

A-T Turner Series

Intended use: To turn patient and to hold patient onto their side.



Product Features

- Fits into the lower back of patient
- Bariatric sizes available
- Available in washable & wipeable fabric
- Latex Free

Options

Sizes: Max Weight / Length

- Offered in both washable (Velcro) or wipeable (coated polyester strap)
 - Regular: 150kg/33"L
 - Bariatric: 300kg/36.6"L

Material

- Polyester Tight Weave Mesh
- Spacer Mesh
- Coated Polyester Webbing - Wipeable
- Polyester Webbing & Velcro - Washable
- Urethane Coated Nylon



Alpha Modalities

A-Pannus/A-PannusSW

Intended use: To assist in moving patient's panniculus away from the body for cleaning skin, wound checks, or dressing changes.



Access port
for cleaning or
exam



Product Features

- Holds panniculus to allow inspection of skin fold areas
- Compatible with loop lift systems or used manually
- SWL 454kg / 1000lbs
- Latex Free

Options

Washable Material-

- Soft inner wicking liner
- Polyester, wicking agent, antimicrobial, polyester webbing

Wipeable Material-

- Ideal for ED and areas that do not have laundry services – wipe with hospital approved disinfectants or 10% bleach solution
- Coated Polyester Webbing, urethane coated nylon, 14oz polyvinyl



A-Limb Series

Intended use: To support limbs for dressing, wound, and/or back care.



Product Features

- Designed to lift and hold limbs eliminating the need to manually hold the limb for dressing changes, pericare, and placing pillows.
- The use of two limb sling will help with patients of size in holding limbs while on their side.
- Available in washable, wipeable, and SPU disposable
- Custom sizes available by request
- SWL 150kg/330lbs
- Latex Free

Options

- A-LimbTWF8: reusable with F8 numerical strap
- A-LimbTWSC: reusable with single clip
- A-LimbSWWL: wipeable with loop or clip
- A-LimbSPU: disposable with loop or clip

Colors:

- Washable available in blue, grey or purple
- Wipeable is only available in blue
- SPU is white



A-SitStandSW

Intended use: Used with a powered mobile lift unit to assist a patient to a standing position.



Product Features

- Addresses issues of having to purchase multiple pieces per size for laundry turn-around.
- Inner safety belt that is adjustable and nonslip.
- Easy to wipe with hospital grade disinfectant and will be ready immediately after cleaning.
- Loop, clip or Arjo-Huntleigh Sara Plus attachment options
- Latex Free

Options

Sizes - color coded tab indicates size

- **Small** (Red) – Inner Safety Belt 29"-46"
- **Medium** (Yellow) – Inner Safety Belt 38"-55"
- **Large** (Green) – Inner Safety Belt 45"-64"
- **X-Large** (Purple) – Inner Safety Belt 55"-73"

Material

- Herculite Fusion III
- 18oz Polyvinyl
- Coated polyester webbing and thread
- Chemical seam seals
- Max Safe Working Load – 300kg



A-WalkVSWWL

Intended use: Ambulation and gait training



Product Features

- Provide safety to patients who are at risk of falling while ambulating
- Multiple sizes identified by color of the vest
- Design allows for toileting by releasing the clip on the groin straps
- Chemical and heat sealed seam
- Max SWL of all vests - 300kg/660lbs
- Latex Free

Options

Sizes-Based upon chest measurement

- **Pediatric** (Blue) 20-26"/50-66cm
- **Small** (Red) 34-38"/86-97cm
- **Medium** (Yellow) 39-43"/99-109cm
- **Large** (Green) 44-48"/112-122cm
- **Xlarge** (Purple) 49-53"/123-135cm
- **Extender Panel** adds 16"/40cm

Material

- Coated polyester webbing and thread
- Metal Cobra buckles
- Nylon coated Urethane
- Spacer mesh polyester



A-BariVest Series

Intended use: Ambulation and gait training



Product Features

- Provide safety to patients of size who are at risk of falling while ambulating
- Multiple sizes identified by color of the vest
- Comes in both washable and wipeable versions
- Cobra Buckles - either plastic or metal
- Scrotum Pouch for comfort
- Latex Free

Options

Sizes - by Max weight

- **Small** (Red) 100kg
- **Medium** (Yellow) 200kg
- **Large** (Green) 300kg
- **Xlarge** (Purple) 400kg


Material

- Polyester
- Coated polyester webbing
- Spacer Mesh
- Urethane Coated Nylon
- Herculite Fusion III



A-WedgeSW

Intended use: Positioning wedge to keep patient on their side or prevent sliding down in bed.

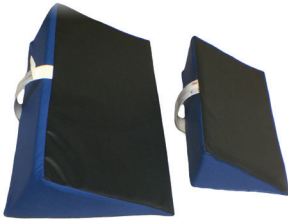

Product Features

- Designed with black non skid surface
- Place between the bed and patient's linen to keep your patient on their side.
- The white wipeable handle allows you to move and position the wedge easily.
- The wedge is sealed to allow for easy cleaning by wiping down with hospital approved disinfectant
- Two size options
- Latex Free

Options

Sizes

- 16"/40.6cm Length
- 23"/58.4cm

Fill Options

- Medical grade foam - firm
- Spun Polyester - soft

Material

- 100% polyester
- Herculite Fusion III
- Waterproof
- Chemical seam sealers



A-BedLadderSW

Intended use: To help patients sit themselves up in bed, assist in moving legs in/out of bed or vehicle and general repositioning while on a flat surface.



Product Features

- MRI Compatible
- Soft, pliable grips
- Can take the place of a trapeze in some cases
- SWL - 272kg/600lbs
- Latex Free

Options

Size

Universal Size 12"/37cm Width 63"/160cm


Material

- Coated polyester webbing and thread
- Herculite Fusion III
- Closed cell foam
- Chemical seam sealer



SPU SERIES

Intended use: Single Patient Use slings that offer the same functions as our TSL, Seat, Limb, Hygiene, and BariVest.

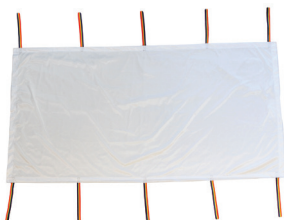

Product Features

- An option for facilities who have challenges with a reusable sling management process
- Variety of sling styles available
- Synthetic breathable material that is strong, comfortable, and can be reprocessed
- Lightweight and individually packaged in clear sleeves
- Can be spot cleaned with hot water if soiled by patient
- **Note:** Slings are packaged 10 per case

Options/Sizes

Flat Repo Sling (TSL)

- 45"X86" / 114x218cm SWL 454kg
- 55"x86" / 140x218cm SWL 500kg

Seated Slings

- M/L Combo SWL 100kg
- X/XX Combo SWL 363kg

Limb Sling

- Universal Size SWL 150kg

Hygiene Sling

- M/L Combo SWL 363kg
- X/XX Combo SWL 363kg


BariKit-SPU

- All-in-one pack includes one of each of the following: Repositioning Sling, Limb Sling, and a BariVest - SWL 454kg

Alpha Modalities

Customizations

Sling Fabric Colors

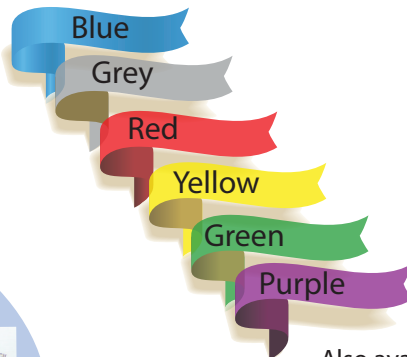

Also available
in white or
other custom
colors



Sling Label Options

Multiple attachment options

Industry
Standard



Color Ladder
Loop

Numerical loop
identification



F8

Aluminum Clip-
Wipe or
wash



Clip

Minimize trip
hazard



Recoil

Wipe or
wash



Mona



Alpha Modalities


Compatibility Statement

All lift motors perform a basic function, which is to go up and down. The interface between the lift, hanger bar and sling is critical for safety and patient comfort.

We test our slings with various manufacturers' hanger bars for compatibility to ensure that they are safe to use and are comfortable for the patient.

We guarantee compatibility with the following manufacturers:

Guldmann AS
Arjo Huntleigh
Prism Medical
Handicare
Tollos
EZ Way
Hill Rom Liko
Etac-Molift

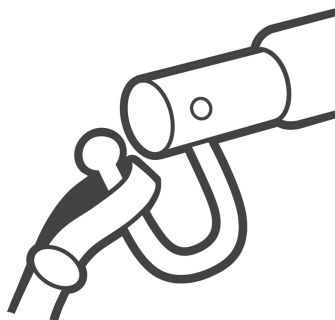

Cleaning instructions for Washable products

- Commercial laundry systems with maximum load 850lbs - tunnel washer voids warranty
- Oxygen based bleach solutions only - chlorine based solution voids warranty
- Fresh water rinse cycle of minimum eight (8) minutes
- Wash temp not to exceed 185F / 85 C
- Dry temp not to exceed 145F / 74C



Cleaning instructions for Wipeable products

- DO NOT LAUNDER - Wipe with hospital approved facility disinfectant for "Non Critical Devices", as per FDA guidelines or 10% bleach solution
- DO NOT soak in solution of chlorine based bleach solution
- DO NOT autoclave or steam clean
- For blood stains: if the blood is fresh and has not dried, *cold water must* be used to clean. DO NOT USE HOT WATER as it sets the stains to iron color from the blood. If it has dried, a spot clean of the stain of 50% Hydrogen Peroxide can be used with cold water.



WARRANTY

Within the first three (3) Years of use of production date, we will repair or replace, at no charge, any sling found to be defective in manufacturing or damaged in use /cleaning. Sling must be returned to us clean, freight prepaid and attached to Return Authorization Number.

ALL PRODUCTS ARE LATEX FREE

