

Your Sling Source for Safe Patient Mobilization

Alpha Modalities



SLING CATALOG

Reusable, Wipeable & SPU Series

September 2019

800.273.5749
www.alphamodalities.com
info@alphamodalities.com



Contents

Guide to Color Coding for Mobility Assessment

Our slings are categorized based on general patient mobility

- RED** **Full Assist:** Patient is unwilling or unable to cooperate
- YELLOW** **Partial Assist:** Patient can follow basic commands but requires assistance
- GREEN** **Stand-by Assist:** Patient can follow basic commands and requires little to no assistance


RED	
Page 1.....	A-TSL Series
Page 2.....	A-Seat Series
Page 3.....	A-SitOnSling
Page 4.....	A-TTurner Series
Page 5.....	A-Pannus Series
Page 6.....	A-Limb Series
YELLOW	
Page 7.....	A-SitStandWE Series
Page 8.....	A-WalkVSSWL
Page 9.....	A-BariV/A-BariVSW Series
GREEN	
Page 10.....	A-WedgeSW Series
Page 11.....	A-BedLadderSW
Page 12.....	SPU Series
Page 13.....	ECO Series
Page 14.....	Customizations
Page 15.....	Compatibility Statement
Page 16.....	Cleaning and Warranty



Alpha Modalities

A-TSL Series

Intended use:
Turns, boosting up in bed, lateral transfers, and cardiac chair transfers.

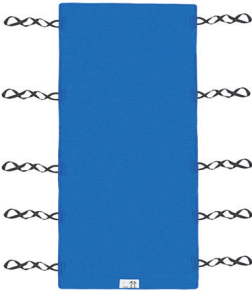

Product Features

- Breathable—designed to remain underneath patient
- Attaches to ceiling and floor lifts w/ 2,4, or 6 point hanger bar systems
- 4 elastic corner mattress straps help prevent sling from moving around on the bed
- Pressure Mapping documentation on low air loss support surfaces
- Available in multiple colors
- Max Load of 454kg /1000lb w/option for 500kg/ 1100lbs

Options

Materials

- Open weave mesh—bathing and general use
- Tight weave mesh—replaces linen on bed
- IT—wicking material w/woven antimicrobial

Sizes

- Width Options: 45", 50", 60", or 66"
- Length Options: 86", 100", or 118"

Attachments

- Color Ladder (CL)—industry standard
- F8 Numerical— numerical identification aids caregivers w/color blindness
- F8 Recoil (RC)- designed to prevent tripping with low support surfaces

Alpha
Modalities

A-SeatSeries

Intended use:

Seated transfers to and from bed, chair, toilet, commode and wheelchair.



Product Features

- Available with or without head support
- Color coded to easily identify size
- Tape sealed seams on wipeable slings
- Max Load 454kg/1000 lbs
- Our Bari Line is made with our IT wicking material and will address patients of size up to 500kg /1100 lbs

Attachments

- Color Ladder (CL)-industry standard
- F8 Numerical- numerical identification aids caregivers w/color blindness
- F8 Recoil (RC)- designed to prevent tripping with low support surface
- Key clip attachment option
- Wipeable Loops- for spray wipeable sling
- Amputee Straps option

Options

Materials

- Open weave mesh-bathing and general use
- Tight weave mesh-replaces linen on bed
- Urethane coated nylon-easy to wipe clean

Sizes/Recommended Weight Ranges:

- **Pediatric** (Dino Print): < 15kg
- **Small** (Red): 14-34kg
- **Medium** (Yellow): 34-66kg
- **Large** (Green): 66-100kg
- **X-Large** (Purple): 100-200kg
- **XX Large** (White): 200-363kg
- **Bari-X** (Grey w purple trim): 363- 500kg

Alpha
Modalities

A-SitOnSling

Intended use:

To provide full body support, without splitting the legs, while transferring to/from bed, chair, toilet, commode and wheelchair.



Product Features

- Our Grey IT material is designed to remain underneath the patient for longer periods
- The edge of sling color is color coded to indicate size
- The sling offers support from the knees to the shoulders on patients
- Placement is accomplished by turning the patient side to side, positioning the sling at the bend of the knee.

Options

Materials

- Grey IT material
- Polyester Webbing and Trim
- Coated Polyester Thread
- Max Load 454kg/1000lbs

Sizes

- **Medium (Yellow Trim)** - Seat width 30in/72cm
- **Large (Green Trim)** - Seat width 32 in/81cm
- **XL (Purple Trim)** - Seat width 35in/89cm
- **XXL (White Trim)** - Seat width 39in/99cm

Alpha
Modalities

A-TTurner Series

Intended use:
To turn patient and hold patient onto their side.



Product Features

- Fits into the lower back of patient
- Bariatric sizes available
- Available in washable & wipeable fabric

Options

Material

- Tight weave mesh
- Spacer mesh
- Coated polyester webbing - Wipeable
- Polyester webbing & Velcro - Washable
- Urethane coated nylon

Sizes: Max Weight / Length

- Offered in both washable (Velcro) or wipeable (coated polyester strap)
 - Regular: 150kg/33"L
 - Bariatric: 300kg/36.6"L

A-Pannus Series

Intended use:

To assist in moving and holding patient's panniculus away from the body for cleaning skin, wound checks, and/or dressing changes.



Access port for cleaning or exam

Optional accessory- O'Ring Headboard Straps



Product Features

- Compatible with loop lift systems or can be used manually
- Comes with an extender belt kit which will extend the length of the bottom attachment loops
- Single Patient Use option available
- Max Load 454kg / 1000lbs

Options

Materials

- IT-wicking material w/woven antimicrobial
- Soft inner wicking liner
- Polyester webbing

Optional Accessory - O'Ring Headboard Straps

- help provide extra lift for patients with a deep panniculus fold and eliminates the need for staff to have to power grip and exert force in order to pull panniculus up and back toward headboard - sold separately

Alpha
Modalities

A-Limb Series

Intended use:

To lift and hold limbs, eliminating the need to manually hold the limb for dressing changes, pericare, and/or placing pillows.



Product Features

- The use of two limbs sling will help with patients of size in holding limbs while on their side.
- Available in washable, wipeable, and SPU disposable
- Custom sizes available by request
- Max Load 150kg/330lbs

Options

Materials

- Tight weave mesh
- Urethane coated nylon

Colors:

- Washable available in blue, grey or purple
- Wipeable is only available in blue
- SPU is white


Attachments

- F8 Numerical- numerical identification aids caregivers w/color blindness
- Key clip attachment option
- Wipeable Loops- for spray wipeable sling
- Rainbow Loop Straps for SPU

A-SitStand Series

Intended use:

Used with a powered mobile lift to assist a patient to a standing position.

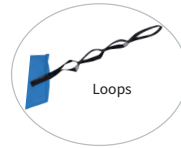

Sara Plus® is a registered trademark of Arjo Inc



Product Features

- Addresses issues of having to purchase multiple pieces per size for laundry turn-around
- Inner safety belt that is adjustable and nonslip
- Easy to wipe with hospital grade disinfectant and will be ready immediately after cleaning
- Max Load 300kg /660lbs

Options

Material

- Outer Belt- Herculite Fusion 9
- Inner Belt- Perfectex
- Wipeable coated polyester webbing on straps
- Coated polyester thread
- Welded seams

Sizes - color coded label indicates size

- Small (Red) – Inner Safety Belt 29"-46"
- Medium (Yellow) – Inner Safety Belt 38"-55"
- Large (Green) – Inner Safety Belt 45"-64"
- X-Large (Purple) – Inner Safety Belt 55"-73"

Attachments

- Loop, clip, or Arjo Sara Plus® options

A-WalkVSSWWL

Intended use:
To assist with ambulation
and gait training.



Product Features

- Provide safety to patients who are at risk of falling while ambulating
- Multiple sizes identified by color of the vest
- Design allows for toileting by releasing the clip on the groin straps
- Chemical and heat sealed seam
- Max Load of all vests - 300kg/660lbs

Options

Material


- Coated polyester webbing and thread
- Metal Cobra buckles
- Urethane coated nylon
- Polyester spacer mesh

Sizes-Based upon chest measurement

- Pediatric (Blue) 20-26"/50-66cm
- Small (Red) 34-38"/86-97cm
- Medium (Yellow) 39-43"/99-109cm
- Large (Green) 44-48"/112-122cm
- Xlarge (Purple) 49-53"/123-135cm
- Extender Panel adds 16"/40cm

A-BariVest Series

Intended use:
To assist with ambulation
and gait training for patients of size.

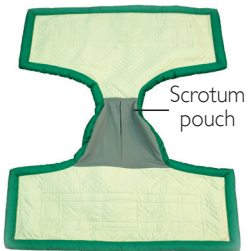

Product Features

- Provide safety to patients of size who are at risk of falling while ambulating
- Multiple sizes identified by color of the vest
- Comes in both washable and wipeable versions
- Cobra Buckles - either plastic or metal
- Washable sling contains a scrotum pouch for comfort
- Now offered in Single Patient Use option

Options

Material

- Polyester
- Coated polyester webbing
- Spacer Mesh
- Urethane coated nylon
- Herculite Fusion III


Sizes - by Max weight

- **Small** (Red) 100kg
- **Medium** (Yellow) 200kg
- **Large** (Green) 300kg
- **Xlarge** (Purple) 400kg

Alpha
Modalities

A-WedgeSW Series

Intended use:
To keep patient on their side or to prevent sliding down in bed.

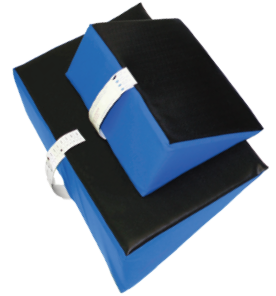

Product Features

- Designed with black non skid surface
- Place between the bed and patient's linen to keep your patient on their side.
- The white wipeable handle allows you to move and position the wedge easily.
- The wedge is sealed to allow for easy cleaning by wiping down with hospital approved disinfectant
- Two size options
- Two fill options

Options

Sizes

- 16"/40.6cm Length
- 23"/58.4cm

Fill Options

- Medical grade foam - firm
- Spun Polyester - soft

Material

- 100% polyester
- Herculite Fusion III
- Waterproof
- Chemical seam sealers

A-BedLadderSW

Intended use:
To help patients sit themselves up in bed, assist in moving legs in/out of bed or vehicle and general repositioning on a flat surface.



Product Features

- MRI Compatible
- Soft, pliable grips
- Can take the place of a trapeze in some cases
- Max Load - 272kg/600lbs

Options

Size
Universal Size 12"/37cm width 63"/160cm

- Material**
- Coated polyester webbing and thread
 - Herculite Fusion III
 - Closed cell foam
 - Chemical seam sealer

A-SPU SERIES

Intended use:
Single Patient Use slings that offer the same functions as our TSL, Seat, Limb, Hygiene, BariVest, and Pannus slings.



Product Features

- An option for facilities who have challenges with a reusable sling management process
- Variety of sling styles available
- Synthetic breathable material that is strong, comfortable, and can be reprocessed
- Lightweight and individually packaged in clear sleeves
- Can be spot cleaned with hot water if soiled by patient
- Slings are packaged 10 per case
- **Now offering a Pannus SPU**

Options/Sizes

Flat Repo Sling (TSL)

- 45"X86" / 114x218cm Max Load 454kg
- 55"X86" / 140x218cm Max Load 500kg

Seated Sling-comes standard w/amputee strap

- M/L Combo Max Load 100kg
- X/XX Combo Max Load 363kg

Limb Sling

- Universal Size Max Load 150kg

Hygiene Sling

- M/L Combo Max Load 363kg
- X/XX Combo Max Load 363kg

BariKit

- All-in-one pack includes one of each of the following: Repositioning Sling, Limb Sling, and a BariVest - Max Load 454kg



Alpha Modalities

A-ECO SERIES

Intended use:
Single Patient Use that can be reprocessed for additional usage.



Product Features

- Synthetic breathable material that is strong and comfortable.
- More cost effective than traditional SPU slings - our quality fabric can be rewashd up to 50 times for additional usage without the worry that the fabric will wear out.
- Variety of sling styles available
- Lightweight and individually packaged in clear sleeves
- Can be spot cleaned with hot water if soiled by patient
- Slings are packaged 10 per case
- 90-day limited warranty

Options/Sizes

Flat Repo Sling (TSL)

- 45'X86" / 114x218cm Max Load 454kg
- 55"x86" / 140x218cm Max Load 500kg

Seated Slings-comes standard w/amputee strap

- M/L Combo Max Load 100kg
- X/XX Combo Max Load 363kg

Limb Sling


- Universal Size Max Load 150kg

Hygiene Sling

- M/L Combo Max Load 363kg
- X/XX Combo Max Load 363kg

Alpha
Modalities

Customizations

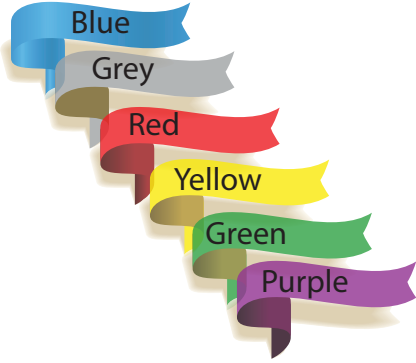


Sling Fabric Colors

Also available in white or other custom colors

Sling Label Options



Multiple attachment options

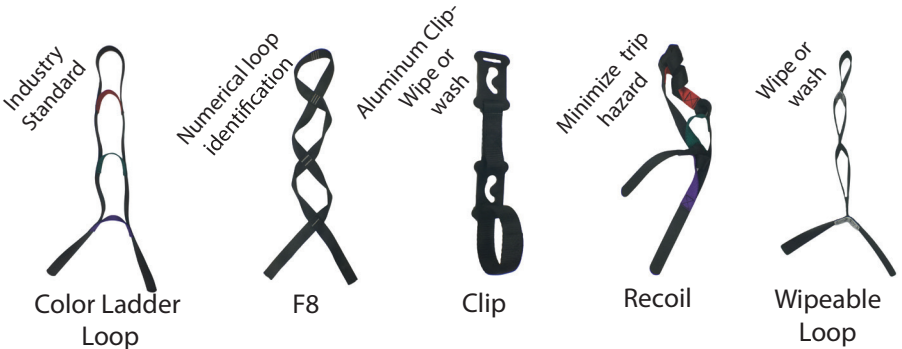

Compatibility Statement

All lift motors perform a basic function, which is to go up and down. The interface between the lift, hanger bar and sling is critical for safety and patient comfort.

We test our slings with various manufacturers' hanger bars for compatibility to ensure that they are safe to use and are comfortable for the patient.

We guarantee compatibility with the following manufacturers:

Guldmann AS
Arjo Huntleigh
Prism Medical
Handicare
Tollos
EZ Way
Hill Rom Liko
Etac-Molift



SLING CLEANING INSTRUCTIONS

Alpha Modalities LLC Washable, Wipeable, SPU and Herculite Fabric products cleaning instructions are listed below. Failure to clean per our guidelines will void product warranty.

Washable products



- Commercial laundry systems with maximum load 850lbs - *tunnel washer voids warranty*
- Oxygen based bleach solutions only - *chlorine based solution voids warranty*
- Fresh water rinse cycle of minimum eight (8) minutes.
- Wash temp not to exceed 185F /85 C
- Dry temp not to exceed 145F / 74C

Wipeable(SW) products



- **DO NOT LAUNDRER** - Wipe with hospital approved facility disinfectant for “Non-Critical Devices”, as per FDA guidelines or 10% bleach solution.
- **DO NOT soak** in solution of chlorine based bleach solution.
- For blood stains: if the blood is fresh and has not dried, **cold water** must be used to clean. **DO NOT USE HOT WATER** as it sets the stains to iron color from the blood. If it has dried, a spot clean of the stain of 50% Hydrogen Peroxide and cold water can be used.

SPU products



- **DO NOT LAUNDRER** - may be spot cleaned with hot water if soiled by patient.

Herculite Fabric Items - (SitStandWE & WedgeSW Series welded products)



- **DO NOT LAUNDRER** - For soils or stains, use neutral soap suds and lukewarm water. **DO NOT use harsh cleansers, solvents or detergents.** For hard to clean spots, use a soft bristle brush and gentle agitation on troublesome spots or stains. Pre-soak heavy, dried-on soil.
- Disinfecting - Quaternary and quaternary/Isopropyl disinfectants are recommended. **DO NOT use phenolic disinfectants.** May be cleaned with a 1:10 dilution of household bleaches containing 5.25% sodium hypochlorite as recommended by the CDC in Atlanta, GA. There is no harmful effect when used in this dilution ratio.
- Disinfectants applied at full concentration or in highly concentrated solutions may damage the fabric
- Iodophor type disinfectants may result in staining.

WARRANTY

Within the first three (3) Years of use of production date, we will repair or replace, at no charge, any sling found to be defective in manufacturing or damaged in use /cleaning. Sling must be returned to us clean, freight prepaid and attached to Return Authorization Number which can be obtained by visiting our Customer Support Center at www.alphamodalities.com.

DO NOT DRY CLEAN OR IRON ANY PRODUCTS
ALL PRODUCTS ARE MADE IN THE USA.
ALL PRODUCTS ARE PROP 65 COMPLIANT AND ARE LATEX FREE.



Specifications and materials are subject to change without notice.

Alpha Modalities
1111 South 344th Street, Suite 205
Federal Way, WA 98003
800.273.5749
www.alphamodalities.com