

Alpha  
Modalities



**Repositioning Slings**  
**Integrated Topsheet (IT) Series**

---


# The Right Sling. The Right Outcome.



Support a healing environment with our proprietary, **moisture wicking and breathable** fabric



Mattress straps may help **minimize pressure points** under the patient

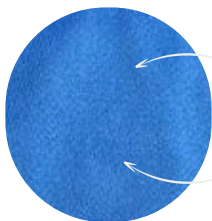

The silver ion embedded fabric is permanent and may help **reduce infection risk**

## Layer Guidance

The IT repositioning sling was designed to replace the fitted sheet. It can remain under patients long term. Follow your local facility guidelines.

### *Best Practice*

**2** layers (IT Sling & Chux)

### *Acceptable Practice*

**3** Layers (Fitted Sheet, IT Sling & Chux)



Contact your Account Manager for testing on breathability, pressure mapping, and radiolucency.

# Customizable Solutions

1

## Choose Your Color

**Green**




**Blue**

**Gray**



2

## Choose Your Sizes

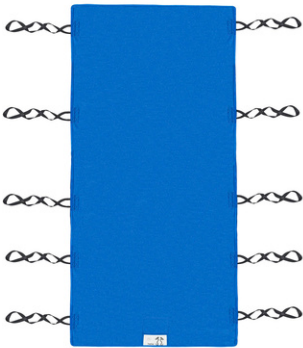

**Blue** or **green** is widely selected for 45" width options.

### Options (inches)

45 X 86  
45 X 100  
45 X 118

50 X 86  
50 X 100  
50 X 118

60 X 86  
60 X 100  
60 X 118

66 X 86  
66 X 100  
66 X 118

**Gray** is recommended to differentiate between smaller sling sizes. Patients of size are most appropriate for 50", 60", and 66" width options.



*\*Additional sizes available upon request.*

# 3

## Choose Your Straps

Color Ladder



Figure 8




Figure 8 Recoil



Figure 8 Recoil was designed with fall prevention and infection prevention in mind. This is the best option for low height patient support surfaces allowing the recoil feature to minimize cross contamination from the floor and reduce risk of trip hazard to caregivers & patients. The numerical identification on both figure 8 options are also beneficial to caregivers with color blindness.

# 4

## Choose Your Label



1111 S. 344th St., Ste 205  
Federal Way, WA 98003  
1-800-273-5749  
www.alpha-modalities.com

**Repositioning & Turning Sling**

Refer to attached Directions for Use label for additional information  
**Maximum/Safe Working Load 1000lbs / 454kg**







**⚠️ This sling must be inspected before each use  
Do Not Use if Damaged - Do Not Repair**

45" / 114cm  
Wide

IDN Logo & Name


IF FOUND, PLEASE CONTACT:

Hospital Address & Contact Information









	Date of First Use MM/YR	
--	-------------------------	--

For more information on use, laundry or maintenance of this sling, go to [www.alpha-modalities.com](http://www.alpha-modalities.com)




Foot End





Label is sewn on the foot end on the side that does not come in contact with the patient.

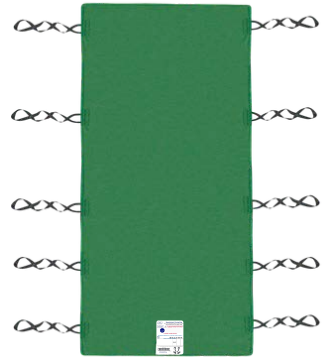

**Prevent Sling Loss**  
with an IDN Specific Label



Available for Mesh and IT Products

*\*Minimum monthly unit spend of 100 units required. Labels are not available on ECO, Wipeable, or Single Patient Use (SPU) Products\**

Example

3

WARRANTY  
YEAR

## Built to last – *Backed for 3 years*

---

### Warranty Coverage

- ✓ Manufacturing defects in materials
  - ✓ Workmanship defects
  - ✓ Loose or unraveled stitching
  - ✓ Structural integrity issues
- ✓ Fabric defects under intended use

When you entrust Alpha Modalities with your safe patient handling and mobility (SPHM) program, we assure you receive an exceptional solution long term. This industry leading warranty period begins either on date of first use or batch tag date, depending upon which is easier to identify. Refer to our website for further details.

***No matter what lift you use,  
we can help.***



We design and manufacture washable, wipeable and SPU slings that are compatible with all ceiling and floor lift hardware sold in North America. With most models, we offer sling options that are compatible with loop and clip sling bars.

# At The Forefront of Sling Innovation

**1** Compatible with all loop style sling bars

**6** Foot end & size indication label with optional IDN customization




**2** Multiple loops on each strap for optimal patient positioning

**3** Mattress straps on every corner to secure sling and reduce bunching

**4** Customizable strap options, sling color, and sizes

**5** Moisture wicking, antimicrobial, breathable IT material


## Safety Starts with Proper Inspection

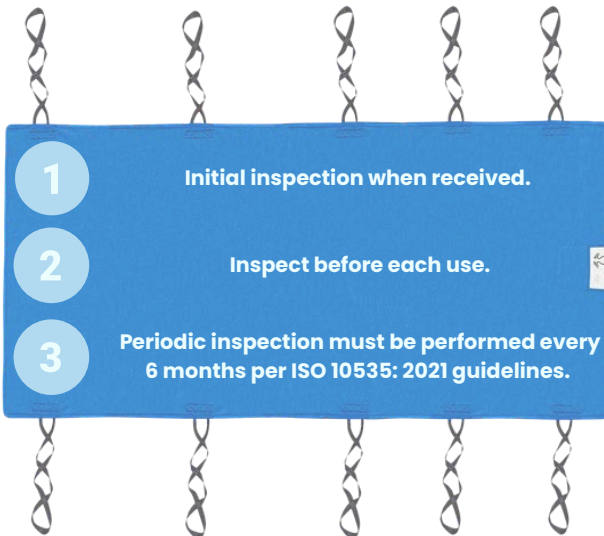

### Immediately remove a sling from service if:

- Any tear, fray, or damage is present
- Stitching is compromised
- Hardware or attachment points are worn or damaged
- Label is missing or unreadable

\*If sling is within warranty period, submit a warranty claim at the Customer Support Center on our website.

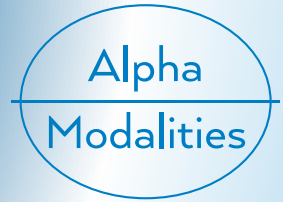


Established in Washington State in 1996, Alpha Modalities has built its reputation on decades of **slings innovation** and manufacturing expertise.



## Our Commitment

Committed to **safety** and **quality**, Alpha Modalities designs slings in alignment with EN ISO 10535 standards. Built for durability, and compatibility across leading lift systems, our solutions are shaped by direct input from caregivers to meet the demands of real-world patient handling.

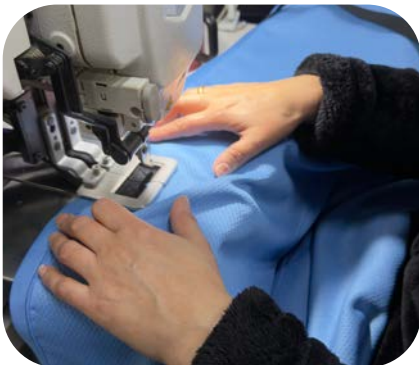

Alpha Modalities' U.S. made, **customizable solutions** reflect the perfect balance of craftsmanship and performance. We look forward to partnering with your facility to advance a sustainable SPHM program.

# Education Support

## Reclined Transfer



## Training Videos



Access our comprehensive training video library

## Turning



## Boost & Supine Transfer



## Quick Reference Guides

Color Ladder Straps  
2 Point Sling Bar

Color Ladder Straps  
4 Point Sling Bar

Figure 8 Straps  
2 Point Sling Bar

Figure 8 Straps  
4 Point Sling Bar




# Technical Specifications

Specifications	TSL-IT Integrated Top (I.T) Sheet
<b>Sizes (inches)</b>	(12) customizable options Width: 45, 50, 60, 66 Length: 86, 100, 118
<b>Washes</b>	185
<b>Material &amp; Color Options</b>	Material: Polyester with antimicrobial and moisture wicking properties (3+) Colors: Blue, Gray, Green, etc.
<b>Safe Working Load (lbs)</b>	1000 - 1100
<b>Warranty</b>	3 years
<b>Mattress Straps</b>	Yes, all 4 corners.
<b>Cleaning</b>	Wash: 185 F Dry:Tumble, Max 145F
<b>Loop Options</b>	(3)Color Ladder, Figure 8, and Figure 8 Recoil; Loop Compatible
<b>Intended Use</b>	2 pt, 4 pt, or 6 pt Loop Style Hanger Bars

## Cleaning Instructions

### IT Series

- Fresh water rinse cycle of minimum eight minutes
- Wash temp not to exceed 185 F / 85 C
- Dry temp not to exceed 145 F / 74C
- Oxygen based solutions only
- \*Chlorine based solutions voids warranty\*
- Commercial laundry systems with Max 850 lbs
- \*Tunnel washer voids warranty\*

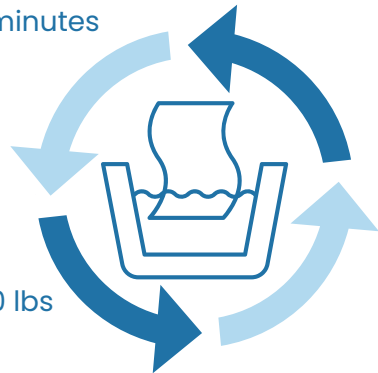

Consult with your Account Manager for sling reprocessing support.

# Can't get enough of the IT? Neither could we!

Explore the full range of IT Solutions at

[www.alphamodalities.com](http://www.alphamodalities.com)



**Pannus**



**Bari Seat**



**Hygiene**



**Sit On**